



**Fastrak®**

**alpha laboratories**  
supplying quality to science

## Fast Racking, Less Packing, Quality Tips

*Fastrak® is probably the most eco-friendly Tip Rack Refill System in the World*

### Safeguard the Quality of Your Work

- Tips Manufactured from Highest Grade Resins
- Excellent Tip-fit Across a Wide Range of Pipettes
- Ultrafine Flexible Tips, Bevelled or Non-bevelled Orifice

**Your Green Choice...Without Compromise**





## The Fastrak Concept...

Fastrak® excels in so many ways...

- **Simple to use** - anyone can do it!
- **Unique 'transfer top'** - ensures contamination is avoided
- **Anti tip-wobble feature** - makes tip mounting easy
- **Small and light** - you get what you need - tips not packaging!
- **Low carbon footprint** - 95% renewable packing material
- **Quality Alpha tips** - made from non-cytotoxic, medical grade, 100% virgin polypropylene - DNA/DNase/RNase free!
- **Unique LaboraTree scheme** - further support for sustainability Fastrak users get trees planted in UK woodlands


## Environmental Features

In all environmental aspects, Fastrak excels. Compared to other popular refill systems, it is lighter, smaller and uses less packaging. In addition, packaging is almost entirely renewable (95%) and far exceeds that of popular brands.




## Transfer Top

Using the unique 'transfer top', you can refill pipette tips without touching the tip or wafer. You can also re-close Fastrak to reduce exposure to contamination on the benchtop.

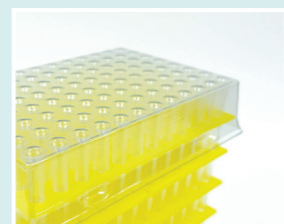

## Anti-wobble Racks

Fastrak tip racks are available with an anti-wobble feature to facilitate tip mounting, especially useful with multichannel pipettes.

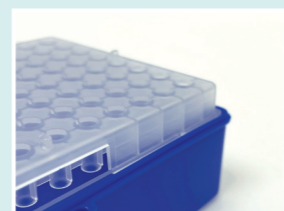

## Minimal Packaging

Fastrak package weight has been minimised to reduce the amount of waste and the fuel needed for transportation.



## The Heart of Fastrak

Transfer top made of recycled PET

Tips are DNA/DNase/RNase free, certified lot by lot

Tip base made from bio-degradable corn derivative

Minimal package weight reduces waste and transport energy consumption

Transfer top keeps tips clean and prevents tip wobble during refilling

You refill without having to touch the tips or racks: so no worries about contamination

When opened Fastrak is extremely stable

Box made from 100% recycled paper

Printed with soya inks to reduce heavy metal pollution

## Rack Features and Benefits

Rack lid opens completely out of the way for pipetting

Patented insert attachment posts automatically lock the refill in place

Robust rack design builds support into the rack, not the insert, so you don't waste plastic from thick inserts

Autoclavable, 96 format racks made from recycled plastic

Hinged lock latch for one-handed operation - it won't pop open if dropped

Rack can be recycled

Exclusive lock hinges won't pop off like side mounted hinges

Every rack lid is labelled with the tip style, lot number, expiry date and volume range

## How Fastrak works

Open Fastrak and grasp the transfer top. Lift straight up.



Align the transfer top over a rack and lower refill into position.



Push down on the transfer top to release the refill plate. Return the transfer top to the Fastrak refill stack and close.



## Quality

A range of features protect the quality of your work



## Green Facts

**Disposal**  
All packaging materials are fully recyclable



Carton



Stack Base



Tip Platform



Transfer Top



# Fastrak® Ordering Information

Fastrak tips are compatible with a wide range of single and multichannel pipettes including:

- Sartorius (Biohit)
- Gilson
- Eppendorf
- Thermo Finnpipette
- Thermo Matrix
- Rainin
- Brand

Naturally it is vital to have a good tip fit for your particular pipette. The ability of a tip to fit a pipette can vary even for the same brand and model due to:

- Tip cone wear
- Manufacturers' use of different tip cones in various geographic locations
- New designs and improvements

The compatibility of our tips with new and existing pipette models is checked on a regular basis and an up-to-date tip fit guide is available on our website.

For Filter Tips with the same pipette compatibility, try our ZAP™ Premier Filter Tip Range. Find out more at [www.alphalabs.co.uk/fastzap](http://www.alphalabs.co.uk/fastzap)

The best way for you to ensure you're completely happy with Fastrak tip fit is to check before you buy.

Please ask us for samples today.









Just call 0800 38 77 32 or email [sales@alphalabs.co.uk](mailto:sales@alphalabs.co.uk)



**Laboratree**

When you join the Laboratree scheme you can make your Fastrak purchase even more green. With each tip pack you will receive a token. Collect tokens and return on the pre-paid card supplied and we'll plant a tree on your behalf. Working in partnership with eforests, we will also supply you with a certificate to mark the occasion. The more tips you use the more trees we plant - it's simple!




		10µl Ultrafine Micro	10µl Ultrafine Extended	200µl Universal	200µl Reference	200µl Ultrafine Reference	300µl Ultrafine	1250µl Ultrafine	20µl Ultrafine LTS*	250µl Ultrafine LTS*	1000µl Ultrafine LTS*
  											
Fastrak®	Non-Sterile	FR1015	FR1010	FR1200	FR2200	FR0250	FR1300	FR1250	FR3020	FR3250	FR3000
	Pack Size	13 x 96	14 x 96	10 x 96	10 x 96	10 x 96	6 x 96	5 x 96	11 x 96	10 x 96	5 x 96
Empty Racks	Square Holes	FR9200	FR9200	FR9200	FR9200	FR9200	FR9000	FR9000	FR9200	FR9200	FR9000
	Round Holes	FR9500	FR9500	FR9500	FR9500	FR9500			FR9500	FR9500	
	Rack Pack Size	10	10	10	10	10	6	6	10	10	6

\*LTS (\*Lite Touch System) is a trademark of Mettler Toledo

**Fast Racking**

**Less Packing**

**Tree Friendly**

JBN205A 05.24